

Awareness Technology, Inc. is proud to introduce our newest fully automated analyzer, ChemWell 2.

The ChemWell 2 is a 2 microwell plate analyzer, performing both ELISA assays and CLIA assays on standard microwell plates or 8 well microwell strips. The flexible analyzer can run mixed modality assays on each of the 2 microwell plates, switching between CLIA and ELISA or dedicating each or both microwell plates to a specific method.

ChemWell 2 is an open system which will allow you to create an assay menu of tests, commercially available from a variety of suppliers. Easy to use software includes, self-checks, patient database, Quality Control and many reporting options.

- Featuring two probes, one for reagents and one for samples, with disposable tip option
- An integrated 8 manifold wash head with adjustable wash volumes, built in aerosol shield
- In-plate mixing with incubation up to 37°C
- Integrated bar code scanning of patient sample IDs
- Sensor-enabled sample rack, which logs and tracks placement of samples that are scanned
- ChemWell 2's powerful software, provides feedback of on-board processes during operation
- Levy-Jennings quality control tracking and reporting


With a proven track record of success, based on reliability, cost effectiveness and innovation, Awareness Technology, Inc. looks forward to working with you.

Designed and built in Florida, USA



ChemWell[®]2



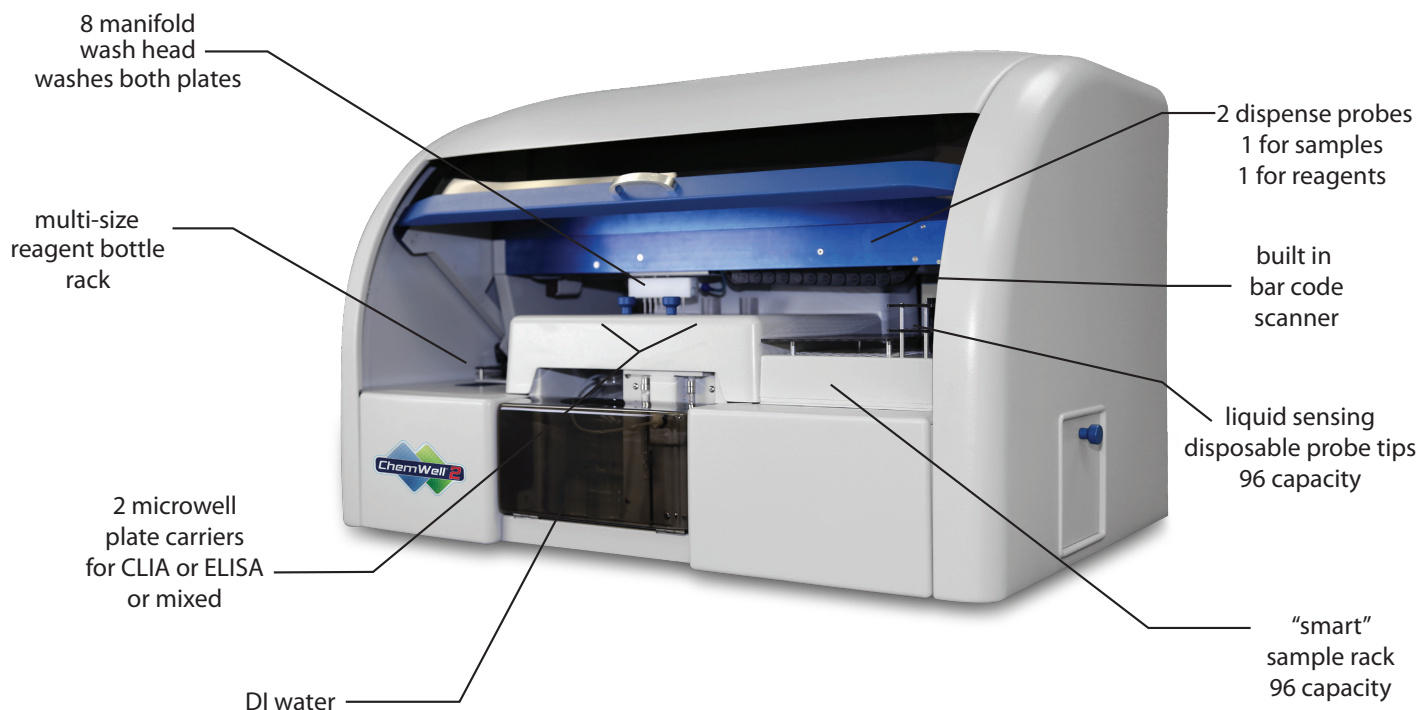

Fluid handling:	Pipetting range 2 μ L – 1.95 mL Typical reaction volume range of 200 μ L to 250 μ L
Probe:	Liquid sensing probe tips monitor remaining volumes, crash protection
CLIA reader:	Linear dynamic range from 1×10^3 to 0.5×10^6 RLU Precision within 4% of target reference over the linear range
ELISA reader:	Accuracy of $\pm 1\% + 0.01$ A or better up to 1.5 A and $\pm 2\%$ or better above 1.5 A
Optical filter range:	4 standard filters 405, 450, 492, 630 nm, optional additional filters available
Standard Sample rack:	96 sample capacity - serum vials or 12x75 mm tubes, "smart" rack detects presence of samples
Standard reagent rack:	33 bottle capacity, mixed sizes
Optional racks:	Reagent racks and sample racks can be manufactured to custom specifications
Bar code scanner:	Integrated bar code scanning with positive patient ID system
Light source:	Halogen lamp
Dimensions:	Height 22 5/8" Width 35 1/8" Depth 23 7/16" Weight ~ 125 lbs



Awareness Technology, Inc. USA & INTERNATIONAL:
1935 S.W. Martin Hwy., Palm City, FL 34990, USA
tel.: +1 772 283 6540, fax: +1 772 283 8020
info@awaretech.com

Awareness Technology, Inc. EUROPE:
Franz-Siegel-Gasse 1, 2380, Perchtoldsdorf, Austria
tel.: +43 (1) 804 81 84, fax: +43 (1) 804 81 85
info@awaretech.eu

Awareness Technology, Inc. MIDDLE EAST:
Q3-146, SAIF Zone, 120543, Sharjah, United Arab Emirates
tel.: +971 (06) 557 80 58, fax: +971 (06) 557 80 59
info@awaretech.com