

sales@awaretech.com

STAT FAX® 4200

Programmable Open System

ELISA Microplate Reader





Stat Fax, serving laboratory needs worldwide for more than 25 years

The Stat Fax® 4200 is a compact, standalone 8-channel microplate reader. Featuring an intuitive touch-screen interface, rapid reading capability of 96 wells, superb optics, six filters with a filter range of 405-700 nm, maintenance-free design.



Features of the Stat Fax® 4200

- ◆ Economical open system & Compact design
- ◆ Test results in 10 seconds per 96 well plate
- ◆ Multiple calculation modes
- ◆ Touch screen display and on-board printer
- ◆ Dependable maintenance-free design
- ◆ Up to 120 user programmable stored tests
- ◆ Custom filter configurations available

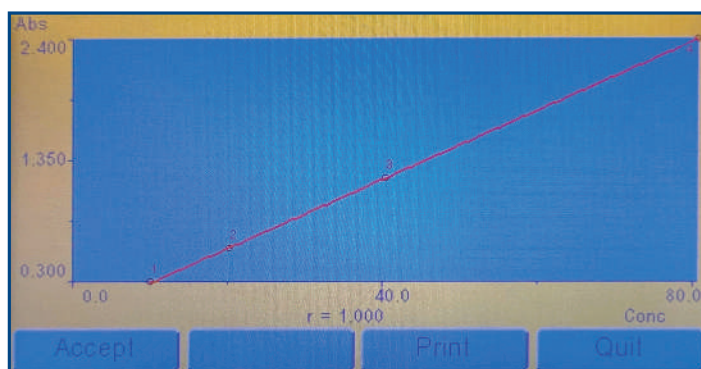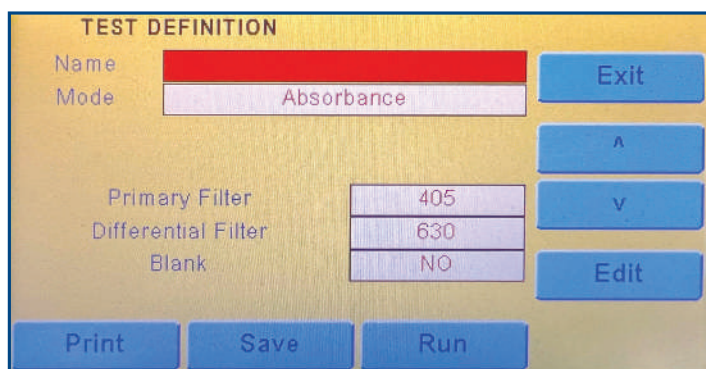

STAT FAX® 4200

Semi-automated ELISA Programmable Plate Reader

Linear Measurement Range	0.0 to 4.0 Absorbance units (A)
Photometric Accuracy	+/- (1% +/- 0.010) 0.0 through 1.5 ABS +/- (2% +/- 0.010) 1.5 through 3.0 ABS
Stability	Drift of no more than 0.005A in 8 hours
Light Source	Halogen lamp with lamp saver feature
Standard Wavelengths	Standard: 405, 450, 492, and 630 nm 6-VIS: 405, 450, 492, 545, 600, and 630 nm Filter options from: 405 - 700 nm
Vessel	Standard 96-well microplates or strip trays
Display	Interactive touch screen 4.5" x 3.5" (14.5 cm) LCD, color graphic display
Printer	Thermal with graphic capability, approximately 64 characters per line
Speed	Reads, calculates and prints results, in approximately 30 seconds per 12-well strip
Calculation Modes	Single point calibration by standard or factor, multi-point calibration with point-to-point curve fit, linear regressions with log, linear, and log-logit selections
Test Storage	Stores 120 tests; stores all parameters including wavelengths, calculations, unit codes, linear and normal ranges, standard values, test names, and previous standard curve.
Dimensions & weight	Approximately 13.5" x 17" x 8.25" (34 cm x 43 cm x 21 cm) weight: 12 lbs (5.5 kg)

ISO 13485:2016



Awareness Technology, Inc. USA & International
tel: +1 772 283 6540 sales@awaretech.com

Awareness Technology, Inc. Europe
tel: +43 (0) 2236 892465 info@awaretech.eu

Awareness Technology, Inc. Asia & Africa
tel: +971 (06) 557 80 58 info@awaretech.com



Dri-Dye® Check Strips

For performance verification of Stat Fax®, ChroMate®,
and other microstrip and microplate readers